

# SECUTEST SIII+ | . . .

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

3-349-397-03  
10/10.11

The SECUTEST SIII+. . . test instruments can be configured for international utilization. The test socket, user interface language and the desired test regulation can all be configured to this end.

**Universal test instrument for testing the electrical safety of:**

- after repair and for periodic testing of electrical devices in accordance with DIN VDE 0701-0702:2008
- after repair and for periodic testing as well as for technical safety checks of electrical medical devices per IEC 62353 (VDE 0751-1) and/or DIN EN 60601 on the basis of the MPG<sup>1)</sup>
- in the production of:
  - electrical measuring, control and laboratory devices per DIN EN 61010
  - data processing systems per DIN EN 60950
  - electrical medical devices per DIN EN 60601
  - household appliances per DIN EN 60335

**DKD Calibration Certificate included**

**Power shutdown as soon as the leakage current exceeds approx. 15mA provides a maximum of safety for the user**

**High-voltage test for:**

**DIN EN 60950, DIN EN 61010, DIN EN 60335 and DIN EN 60601**



Device with optional SI module

## Features

**The device under test can be connected:**

- to the test socket with or without adapter for various types of mains connection
- to the connector jacks if the device under test does not have a mains plug
- with an adapter for extension cables with or without multiple outlet sockets
- connection for BE<sup>2)</sup> and FE<sup>3)</sup>
- 10 application parts can be connected individually or in groups

**Mains Plug Polarity Reversal**

Mains plug polarity need not be reversed manually. Polarity reversal is accomplished internally during the test sequence (except if 3-phase current adapters are used).

**Automatic Recognition**

of mains connection errors and protection class (I or II). Measurement is automatically disabled in the event of danger.

**Display**

Menus, setting options, measurement results, instructions and error messages, as well as online help and schematic diagrams for test setups, can all be displayed at the backlit, dot matrix LCD.

<sup>1)</sup> MPG = German medical product law  
<sup>2)</sup> BE = Operational earth  
<sup>3)</sup> FE = Functional earth

**Menu Driven Test Sequences**

Fully automatic or manual

**Selectable Test Current for Protective Conductor Testing (4-pole Measurement)**

With 200 mA, 10 A or 25 A test current

**Insulation Test**

By means of insulation resistance or equivalent leakage current measurement or high-voltage test

**Leakage Current Test**

With measurement of earth, housing or patient leakage current, patient auxiliary current, protective conductor current, contact current, verification of absence of voltage by means of current measurement or device leakage current

**Basic Instrument and Expansion Features**

The test instrument can be configured with specific features required for the given application (see table on page 6).

**Data Interface for PC, Printer and Barcode**

**Expandable**

The SECUTEST SI+ option expands the basic instrument into a unique data logger with memory and alphanumeric keypad for data entry. All required reports can be generated, and data can be analyzed and managed with the help of user-friendly WINDOWS software.

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

### Applications

#### Testing for the Electrical Safety of Electrical Equipment in Accordance with BGV A3

The test instrument can be utilized for quick and safe testing of repaired or modified electrical devices and as well as for periodic testing in accordance with DIN VDE 0701-0702:2008.

The following are measured in accordance with the standards:

- Protective conductor resistance
- Insulation resistance
- Protective conductor current for SC1 devices
- Contact current for SC2 devices
- Absence of voltage at exposed conductive parts (= contact current)

Measuring methods for leakage current measurements:

- Direct measurement
- Equivalent leakage current
- Residual current

#### Testing for the Electrical Safety of Electrical Medical Devices in Accordance with the German Medical Product Law (MPG) and the associated Operator's Regulations

The test instrument with feature KA01 is used for quick and safe testing and measurement of repaired or modified electrical medical devices or their components (e.g. patient ports) in accordance with IEC 62353/DIN EN 62353 (VDE 0751) and EN 60601, see also **SECULIFEST** oder **SECULIFE ST HV**.

Observance of technical safety requirements allows the user of the test instrument to operate electrical medical devices in a hazard-free fashion. The safety of the patient is also assured through the use of tested electrical medical devices.

#### The following are measured in accordance with IEC 62353 (VDE 0751-1) regulations:

- Protective conductor resistance
- Insulation resistance
- Equivalent device leakage current
- Equivalent patient leakage current
- Device leakage current
- Patient leakage current (AC/DC portions are measured separately)

Measuring methods for leakage current measurements:

- Direct measurement
- Equivalent leakage current
- Residual current

#### A software upgrade (optional) allows for measurement in accordance with EN 60601 regulations, see features on page 6 (with the following single-fault conditions: voltage at application part, interrupted neutral and interrupted protective conductor, with automatic polarity reversal L-N)

- Protective conductor resistance
- Insulation resistance
  - L and N connected to protective conductor
  - Application parts connected to protective conductor
- Earth leakage current, housing leakage current, patient leakage current, patient auxiliary current

#### The following additional test conditions can be selected:

- Interrupted operational earth  $\perp$  and equipotential bonding  $\nabla$
- Housing to ground, application parts to ground

#### Function Test with Power Analysis

(also suitable for high power devices under test up to 16 A)

The device under test can be subjected to a function test with mains voltage via the integrated test socket.

The following are measured or automatically calculated during the function test:

- Line voltage
- Residual current
- Power consumption
- Active and apparent power
- Power factor
- Electrical energy
- On-time

#### Multimeter Functions

Extensive multimeter functions including temperature measurement expand measuring options for the user in a sensible fashion. The following individual measurements can be performed:

- Direct and alternating voltage (momentary and min/max values)
- Resistance
- Voltage against PE, e.g. phase detection
- Current and protective conductor resistance with clip-on meter (accessory)
- Temperature with Pt100 or Pt1000 (accessory)

#### High-Voltage Test with Direct Voltage (SECUTEST SIII+ H)

The mains plug of the device under test (safety class I and II devices) is connected to the test socket at the test instrument.

The test instrument monitors the mains connection. Incorrect or dangerous mains connection is indicated, and measurement is disabled in the event of danger.

Use of the test instrument for high-voltage testing is trouble-free because DIN VDE 0104 does not apply. The high-voltage test is performed with direct voltage. In order to comply with requirements for alternating voltage, testing is performed with 1.5-fold direct voltage. This multiplying factor is applied automatically during testing.

This DC high-voltage test complies with EN 60601 3rd edition/ EN 50106 (VDE 0700 part 500), as well as with other standards.

#### Report Functions

All values required for electrical device approval reports or device log books (e.g. for ZVEH) can be measured with the test instrument.

All measured data can be documented and archived with the measurement and test report, which can be saved to memory and printed out from a PC.

The measurement and test report substantiates regular maintenance and testing for users of electrical devices.

The SECUTEST SI+ module (accessory equipment), a memory with integrated interface and keypad which can be mounted inside the lid of the test instrument, expands the applications range of the test instrument.

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

The test instrument has been manufactured and tested in accordance with the following standards:

IEC/EN 61 010-1:2001 VDE 0411-1:2002	Safety requirements for electrical measurement, control and laboratory devices – General requirements
DIN VDE 0404 Part 1: 2002	Test and measuring equipment for testing the safety of electrical devices – General requirements
DIN VDE 0404 Part 2: 2002	– Testing equipment for tests after repair, modification or in the case of periodical tests
DIN VDE 0404 Part 3: 2005	– Equipment for periodical tests and tests prior to commissioning medical electrical devices or systems
DIN EN 60 529/ VDE 0470 Part 1	Test instruments and test procedures, protection provided by enclosures (IP code)
DIN EN 61 326-1 VDE 0843-20-1	Electrical equipment for measurement, control and laboratory use – EMC requirements – Part 1: General requirements

### Regulations for the Use of the Test Instruments

Devices under test to be tested in accordance with the following regulations	Testing after Repairs / Periodic Testing		Routine Testing			
	DIN VDE 0701-0702:2008	IEC 62353:2007 DIN EN 62353:2008 (VDE 0751-1)	DIN EN 60950	DIN EN 61010	DIN EN 60335/EN 50106	IEC 60601/DIN EN 60601 *
Electrical devices	•			•		
Appliances and electric equipment	•				•	
Mains operated electronic devices	•					
Hand-held electric tools	•					
Extension cables	•					
Data processing devices	•		•			
Electrical medical devices, application parts		•				•

\* only test instruments with Feature KA01


MPG	German medical product law: Regulation for the setup and operation of active medical products
-----	--

### Table: Individual Measurements and Regulations

Individual Measurements per Regulation	Test Current [A]	DIN VDE 0701-0702	DIN VDE 0701 Part 1	DIN VDE 0701 Part 2/40	DIN VDE 0701 attachment E	DIN EN 60950	DIN EN 61010	DIN EN 60335	IEC 62353 (VDE 0751-1)	IEC 601/EN 60 601 2nd	IEC 601/EN 60 601 3rd
		Protective Conductor Resistance	0.2 10 25	•	•	•					•
Insulation Resistance		•	•								
Equivalent Leakage Current		•	•	•							
High-Voltage Test						•	•	•		AC	AC
Equivalent (Device) Leakage Current								•	•		
Equivalent Patient Leakage Current									•		
Residual Current		•	•						•		
Contact Current		•	•								
Absence of Voltage (exposed conductive parts)		•		•							
Housing Leakage Current						•	•			•	•
Earth Leakage Current										•	•
Patient Leakage Current									•	•	•
Total Patient Leakage Current										•	•
Patient Auxiliary Current										•	•
Device Leakage Current									•		
Single Fault Conditions	N PE Mains at Application Part						•	•		•	•

Sample displays, menu-driven operation:


#### Test Sequence Setup

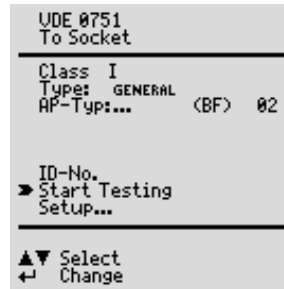

#### Visual Inspection



#### Protective Conductor Test

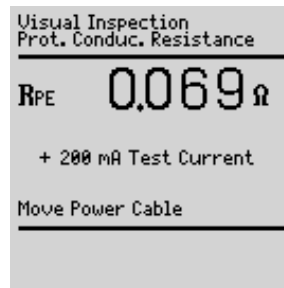


#### Leakage Current Measurement




Sample displays, online help:


#### Online Help Texts

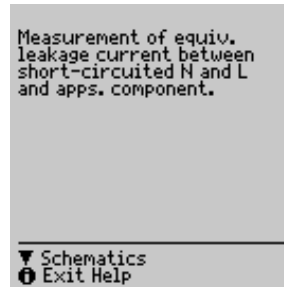


#### Technical Data

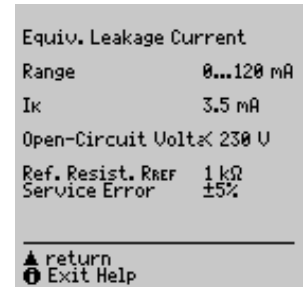

#### Schematic Diagrams

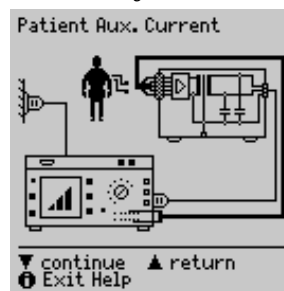


#### Internal Circuit Diagrams

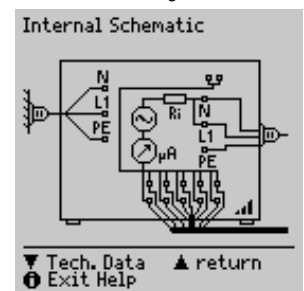


Sample reports with measurement results:

#### Test Sequence Results

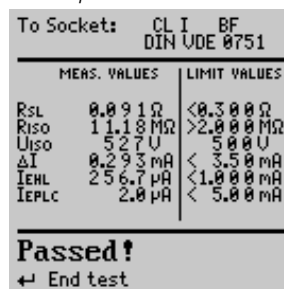

#### Function Test Results

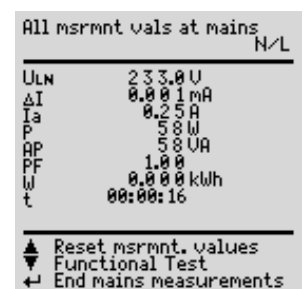

### Characteristic Values

	Measured Quantity	Measuring Range/ Nominal Range of Use	Resolu- tion	Nominal Voltage $U_N$	Open- Circuit Voltage $U_0$	Nominal Current $I_N$	Short- Circuit Current $I_K$	Internal Resis- tance $R_i$	Refer- ence Resis- tance $R_{REF}$	Measuring Uncertainty <sup>3)</sup>	Intrinsic Uncertainty <sup>3)</sup>	Overload Capacity	
												Value	Time
Refer to page 3 for assigning individual measurements to the regulations	Device Protective Conductor Resistance $R_{PE}$	0.000 ... 2.100 $\Omega$	1 m $\Omega$	—	4.5 ... 9 V DC	—	> 200 mA DC	—	—	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$ > 10 d	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont.
		2.11 ... 31.00 $\Omega$	10 m $\Omega$	—	< 6 V AC	—	> 10 A AC <sup>4)</sup> > 5 s	—	—			no protection <sup>5)</sup>	
	Insulation Resistance $R_{ISO}$	0.050 ... 1.500 M $\Omega$	1 k $\Omega$	50 ... 500 V DC	1.0 • $U_N$ ... 1.5 • $U_N$	> 1 mA	< 10 mA	—	—	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont.
		1.01 ... 10.00 M $\Omega$	10 k $\Omega$										
	Equivalent Leakage Current $I_{EL}$	0.00 ... 21.00 mA	10 $\mu$ A	—	230 V ~ -20/ +10 %	—	< 3.5 mA	> 72 k $\Omega$	2 k $\Omega$	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont.
		20.1 ... 120.0 mA	100 $\mu$ A										
	Contact Current (Absence of Voltage) $I_{probe}$	0 ... 3.500 mA	1 $\mu$ A	—	—	—	—	2 k $\Omega$	—	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont. <sup>2)</sup>
	Residual Current $I_{DI}$ between L and N	0.000 ... 3.100 mA ~	1 $\mu$ A	—	—	—	—	—	—	$\pm(10\% \text{ rdg.} + 10 \text{ digits})$ > 10 digits	$\pm(5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	1)	1)
		3.00 ... 31.00 mA ~	10 $\mu$ A										
	Equivalent Device and/or Equivalent Patient Leakage Current $I_{EDL}$ and/or $I_{EPL}$	0.0 ... 310.0 $\mu$ A	0.1 $\mu$ A	—	230 V ~ -20/ +10 %	—	< 3.5 mA	> 72 k $\Omega$	1 k $\Omega$ $\pm 50 \Omega$	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont. <sup>1) 3)</sup>
		0.000 ... 2.100 mA	1 $\mu$ A										
		2.101 ... 21.00 mA	10 $\mu$ A										
	Leakage Current $I_L$ <sup>2)</sup>	0.0 ... 310.0 $\mu$ A	100 nA	approx. line voltage <sup>6)</sup>	—	—	—	1 k $\Omega$	—	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$	$\pm(2.5\% \text{ rdg.} + 5 \text{ digit})$ > 10 digit	253 V	cont. <sup>1) 3)</sup>
All Leakage <sup>7)</sup> Current $I_L$		0.210 ... 3.600 mA	1 $\mu$ A										
		3.10 ... > 15.00 mA	10 $\mu$ A										

Function	Measured Quantity	Measuring Range / Nominal Range of Use	Resolu- tion	Open- Circuit Voltage $U_0$	Short- Circuit Current $I_K$	Internal Resis- tance $R_i$	Measuring Uncertainty	Intrinsic Uncertainty	Overload Value	Capacity Duration		
Functions Test	Nominal Voltage $U_{L-N}$	103.5 V ... 126.5 V 207.0 ... 253.0 V ~	0.1 V	—	—	—	—	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$	253 V	cont.		
	Load Current $I_V$	0 ... 16.00 A RMS	10 mA	—	—	—	—	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$	20 A	10 min		
	Active Power P	0 ... 3700 W <sup>9)</sup>	1 W	—	—	—	—	$\pm(5\% \text{ rdg.} + 10 \text{ digits})$ > 20 digits	253 V	cont.		
	Apparent Power S	0 ... 4000 VA	1 VA	Calculated Value $U_{L-N} \cdot I_V$				$\pm(5\% \text{ rdg.} + 10 \text{ digits})$ > 20 digits	20 A	10 min		
	Power Factor PF, sinusoidal: $\cos \varphi$	0.00 ... 1.00	0.01	Calculated Value P / S, Display > 10 W				$\pm(10\% \text{ rdg.} + 5 \text{ digits})$				
	Residual Current $\Delta I$ between L and N	0.00 ... 31.00 mA ~	10 $\mu$ A	—	—	—	—	$\pm(10\% \text{ rdg.} + 10 \text{ d})$ > 10 digits	$\pm(5\% \text{ rdg.} + 5 \text{ digits})$	1)	1)	
$U_{AC/DC}$	Voltage	0 ... 253.0 V —, ~ and —	0.1 V	—	—	—	—	$\pm(5\% \text{ rdg.} + 10 \text{ d})$	253 V	cont.		
$U_{probe}$	Probe Voltage	0 ... 253.0 V —, ~ and —	0.1 V	—	—	—	—	$\pm(2.5\% \text{ rdg.} + 5 \text{ digits})$ > 10 digits	253 V	cont.		
R	Resistance	0 ... 150.0 k $\Omega$	100 $\Omega$	< 20 V -	1.1 mA	—	—	$\pm(1\% \text{ rdg.} + 3 \text{ digits})$	253 V	cont.		
$I_{clip}$	Current via Clip-On Current-Voltage Converter WZ12C	0.000 ... 10.00 A ~	1 mA	—	—	1.5 M $\Omega$	—	$\pm(3\% \text{ rdg.} + 10 \text{ digits})$ > 10 digits	253 V	cont.		
		0 ... 100 A ~	1 A	—	—	1.5 M $\Omega$	—	without clip	253 V	cont.		
Temp	Temperature with Pt100 / Pt1000 Sensor	-200 ... -50 °C	1 °C	< 20 V -	1.1 mA	—	—	—	—	$\pm(2\% \text{ rdg.} + 1 \text{ } ^\circ\text{C})$	10 V	cont.
		-50.1 ... +300.0 °C	0.1 °C							$\pm(1\% \text{ rdg.} + 1 \text{ } ^\circ\text{C})$	10 V	cont.
		+300 ... +850 °C	1 °C							$\pm(2\% \text{ rdg.} + 1 \text{ } ^\circ\text{C})$	10 V	cont.

- 1) As of 25 mA: shutdown by residual current measurement within 100 ms
- 2) Except for earth leakage current: only 0.000 to 3.100 mA
- 3) Measuring circuit is highly resistive, indication at display
- 4) Measurement with AC test current is not possible at jacks (1) to (3); feature G01: > 25 A: Short-circuit current is less than 25 A if the SK5 special cable is used.
- 5) Test duration max. 40 s, protection against overheating: measurement cannot be restarted until a waiting period of 1 minute has elapsed.
- 6) Calculated value: max. 253 V
- 7) AC and DC are measured for patient leakage current and patient auxiliary current.
- 8) The data are only valid for the values displayed at the test instrument. Data which are transmitted via the RS232 interface may deviate.
- 9) Measured value P and calculated value S are compared, and the smaller value is displayed.

Key: rdg. = reading, d = digits  
 $I_L$  = patient, housing and earth leakage current, as well as patient auxiliary current

### High Voltage Test (feature F02 or SECUTEST SIII+ H)

#### Transducer

Nominal Voltage, AC	$U_N$ - adjustable in 10 V steps in 100 V steps	0.5 ... 0.99 kV 1 ... 4 kV $((U_N \cdot 1.5) \cdot 1.011) + 60 \text{ V}$
Open-Circuit Voltage, DC	$U_0$	$\pm 1.5\%$
Intrinsic Uncertainty, $U_0$	$U_0$	$\pm 1.5\%$
Nominal Current	per DIN VDE 0104	< 3.5 mA DC
Short-Circuit Current	discharge current from 6 x 2.7 nF	> 5 A at 5 kV
Resistance to Interference Voltage		none

#### Measuring

Measuring Range	Display Range	Intrinsic Uncertainty, $U_0$
0 ... $U_{omax}$	0.000 ... > 10.00 kV DC	$\pm 1.5\% \text{ rdg.} + 2 \text{ digits}$

# Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

## Testing for Correct Mains Connection

The test instrument automatically recognizes mains connection errors, if the conditions in the following table have been fulfilled. The user is informed of the type of error, and all measuring functions are disabled in the event of danger.

Type of Mains Connection Error	Message	Condition	Measurements
Voltage at protective conductor PE to finger contact	Text appears at LCD	Press  key $U > 40 \text{ V}$	disabled
Protective conductor PE and phase conductor L reversed and/or neutral conductor N interrupted	lamp lights up	Voltage at PE $> 65 \text{ V}$	impossible (no supply power)
Contact voltage at protective conductor PE to neutral conductor N or phase conductor L	Text appears at LCD	$U > 25 \text{ V}$	disabled, although disabling can be deactivated (e.g. IT network)
Mains voltage too low	lamp lights up	$U_{L-N} < 90/180 \text{ V}$	possible under certain circumstances

## Influencing Quantities and Influence Error

Influencing Quantity/ Sphere of Influence	Designation per DIN VDE 0404	Influence Error $\pm \dots \%$ of Measured Value
Position Change	E1	—
Change in Test Setup Supply Power	E2	2.5
Temperature Fluctuation	E3	Specified influence error applies per 10 K change in temperature: 1 in case of PE measurement 0.5 of all other measuring ranges
0 ... 21 °C and 25 ... 40 °C		
Current at Device Under Test	E4	2.5
Low-Frequency Magnetic Fields	E5	2.5
Impedance at Device Under Test	E6	2.5
Capacitance, Insulation Measurement	E7	2.5
Waveshape of Measured Current	E8	2 for capacitive load (for equivalent leakage current)
49 ... 51 Hz		1 (for contact current)
45 ... 100 Hz		2.5 for all other measuring ranges

## Reference Ranges

Line Voltage	115/230 V $\pm$ 0.2 %
Line Frequency	50/60 Hz $\pm$ 0.1 %
Waveshape	sine (deviation between effective and rectified value < 0.5%)
Ambient Temperature	+23 °C $\pm$ 2 K
Relative Humidity	40 % ... 60 %
Load Impedance	linear

## Nominal Ranges of Use

Line Voltage	103.5 V ... 126.5 V or 207 V ... 253 V
Line Frequency	50 Hz or 60 Hz
Line Voltage	
Waveshape	sine
Temperature	0 °C ... + 50 °C

## Ambient Conditions

Storage Temperature	- 20 °C ... + 60 °C
Operating Temp.	- 10 °C ... + 50 °C
Accuracy Range	0 °C ... + 50 °C
Relative Humidity	max. 75%, no condensation allowed
Elevation	max. 2000 m

## Electromagnetic Compatibility

Product standard DIN EN 61326-1

Interference emission		Class
EN 55011		B
Interference immunity	Test Value	Evaluation Criteria
EN 61000-4-2	Contact/Atmos. - 4 kV/8 kV	A
EN 61000-4-3	3 V/m or 1 V/m	A
EN 61000-4-4	1 kV	B
EN 61000-4-5	1 kV bzw. 2 kV	A
EN 61000-4-6	3 V/m	A
EN 61000-4-11	0.5/1/25 Periods	A
	250 Periods	C

## Power Supply

Line Voltage	103.5 V ... 126.5 V or 207 V ... 253 V
Line Frequency	50 Hz or 60 Hz
Power Consumption for 10 A test current	approx. 30 VA
for 25 A test current for function test	approx. 95 VA, test duration max. 40 s approx. 180 VA, test duration max. 40 s continuous max. 3600 VA, power is conducted through the instrument only, switching capacity $\leq$ 16 A

## RS 232 Data Interface

Type	RS 232C, serial, per DIN 19241
Format	9600, N, 8, 1
Connector	9-pin subminiature socket connector

## Electrical Safety

Safety Class	I per IEC 61010-1/EN 61010-1/ VDE 0411-1
Nominal Voltage	115/230 V
Test Voltage	3.7 kV 50 Hz
Measuring Category	250 V CAT II
Pollution degree	2
Safety Shutdown	for residual current at device under test > 25 mA, disconnecting time < 100 ms probe current > 10 mA, < 1 ms

## Mechanical Design

Display	multiple backlit dot matrix display, 128 x 128 pixels
Dimensions	test instruments without high-voltage module: LxWxH: 292 mm x 138 mm x 243 mm test instruments with high-voltage module: LxWxH: 292 mm x 138 mm x 300 mm
Weight	standard device: approx. 4.5 kg device with HV test: approx. 5.24 kg device with 25 A PE test: approx. 5.5 kg with 25 A PE and HV test: approx. 5.9 kg
Protection	housing: IP 40, connections: IP 20 per DIN VDE 0470 Part 1/EN 60529, Extract from table on the meaning of IP codes

IP XY (1 <sup>st</sup> digit X)	Protection against foreign object entry	IP XY (2 <sup>nd</sup> digit Y)	Protection against the penetration of water
2	$\geq 12.5 \text{ mm } \varnothing$	0	not protected
4	$\geq 1.0 \text{ mm } \varnothing$	0	not protected

# SECUTEST SIII+ | . . .

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

### Standard equipment for SECUTEST SIII+. . .

- |   |                                   |
|---|-----------------------------------|
| 1 test instrument   | 1 Calibration Certificate per DKD |
| 1 probe cable with test probe depending upon test instrument features | 1 operating instructions          |
| 1 plug-on alligator clip for test probes                              | 1 carrying strap                  |
| 3 plug-on quick-connect terminals                                     |                                   |

### Features and Options

List of possible options:

Feature		00	01	02	03	04	05	06	07	08	09	10	11	XX
Mains Connection for Country of Use	B	D	D + ser-vice socket	UK <sup>4)</sup>	F/CZE		DK <sup>4)</sup>			China/AUS <sup>4)</sup>	CH		Adapter kit <sup>2) 4)</sup>	
User Interface Language	C	D	UK	F	I	E	CZE	NL						
High-Voltage Test HV DC	F	without		max. 6,126 kV DC (≤ 4 kV AC)										
AC Test Current 50/60 Hz for Protective Conductor Measurement	G	10 A	25 A	without										
10 + 2 Sockets for Application Parts	J	without	with											
Test Sequence for IEC 60601	KA	without	with <sup>3)</sup>											
Data Memory for up to 125 Tests <sup>5)</sup>	KB	without	with											
Recognition of Probe on Protective Conductor	KD	without	with											
Direct Printing after each Measurement for Auto-matic Test Sequences <sup>1)</sup> via RS232	KE	without	with											

<sup>1)</sup> Each measured value is documented in this case, as opposed to the results of a test sequence for which the poorest value for each given test is displayed (via the PSI module, the SECUSTORE memory adapter or a PC)

<sup>2)</sup> Adapter kit for international use (equipped with Feature B01)

<sup>3)</sup> only possible with Feature J01

<sup>4)</sup> for mains connections B02, B05, B08 and/or if adapter (feature B11) is applied: HV-DC max. 1.5 kV DC

<sup>5)</sup> without function test values and without comments on DUT

Enter the designation of the basic instrument to your order, i.e. M7010, and only those desired features which are other than 00!

Example of a complete type designation (= article number, = order designation) for a SECUTEST SIII+. . . :

SECUTEST SIII+. . . with Swiss plug and Swiss socket in French language, without high-voltage test, with AC test current 25 A, without sockets for application parts, without test sequence for IEC 60601, without data memory, with recognition of probe on protective conductor, without direct print-out option:

Features: **M7010 B09 C02 F00 G01 J00 KA00 KB00 KD01 KE00**

### Standard Types

Type	Designation / Feature Combination	Article number
SECUTEST SIII+ H	including selectable test current ±200 mA DC or <b>25 A AC (G01)</b> , including <b>high voltage testing up to 6 kV DC (F02)</b> including sequences for IEC 61010, IEC 60335, IEC 60950, including data memory for up to 125 tests <sup>5)</sup> (KB01)	M7010-V013

Features which are additionally required can be retrofitted by GMC-I Service GmbH upon request.

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

### Feature KA01: Tests per IEC 60 601/EN 60601

Measurements in accordance with this standard are made possible by uploading the appropriate software to the instrument with the help of a PC via the included interface cable. Special features:

- Patient ports can be assigned to groups
- Automatic sequence under all single-fault conditions

### Feature KB01 (in Standard Model SECUTEST SIII+ H):

#### Data memory for up to 125 tests Memory expansion for test results

If no (P)SI module is connected, up to **125 test results\*** are stored to the test instrument. The test results can be viewed on the display and printed out, for example, via a terminal program. The test results are sorted in chronological order and are shown with the associated ID number. If no ID number has been allocated, date and time of storage are automatically saved instead. Alternatively, a consecutive number can be entered.

\* without function test values and without comments on DUT

#### Storage of parameter settings for test sequences per IEC/DIN EN 60335/60950/61 010

Test sequences can be configured on-site and performed in the appropriate selector switch position in accordance with the respective regulations. These **configurations for various test sequences** are stored to the test instrument and can be reactivated as required.

### Feature KD01: (Recognition of Probe on Protective Conductor) (in Standard Model SECUTEST SIII+ H)

For freely configurable instruments the feature additionally includes a 5 m long probe cable with a test probe. The protective conductor measurement is expanded to include the function: "automatic recognition of measuring point change".

During protective conductor measurement, the instrument recognizes whether or not the probe is in contact with the protective conductor, and indicates these two possible conditions by means of acoustic signals. This function is helpful if several protective conductor connections need to be tested.

### Feature KE01: Direct Print-Out (in Standard Model SECUTEST SIII+ H)

After completion of each test (individual test or at the end of a test sequence), test results are read out directly via the RS232 interface.

**User interface languages** which are not included as a standard feature can also be uploaded from our homepage ([www.gossen-metrawatt.com](http://www.gossen-metrawatt.com)). One language at a time can be uploaded to the test instrument.

### Software Installation Requirements

#### Software:

– MS WINDOWS 2000 or XP.

#### Hardware:

- IBM compatible WINDOWS PC, as from 200 MHz Pentium-CPU with at least 64 MB RAM
- SVGA monitor
- Hard disk with at least 20 MB available memory
- MICROSOFT compatible mouse

## Accessories

### Memory and Input Module SECUTEST SI+

Values measured by the test instrument can be stored to this module, and can be furnished with comments with the help of the alphanumeric keypad. The LCD panel at the test instrument is used as a display for the module. Statistical analysis of the measurement results is also possible (percentage of tests which have been successfully passed). The SI module is screwed into the lid of the test instrument in a space-saving fashion.



Please request our SECUTEST SI+ data sheet for further information.

### SECUSTORE – Memory Adapter for SECUTEST...

Test reports and individual test steps and/or test series can be directly written („printed“) from the SECUTEST... test instrument to the memory adapter. They can be subsequently read out and processed at the PC. This memory adapter is therefore an excellent tool for the archiving and transmission of test reports.

Depending upon the size of the reports, up to 1,000 test reports or test series can be stored.

Please note that the direct print-out option must be enabled in the SECUTEST... to allow for the storage of individual test steps and test series.

### Comparison of Memory Adapters / Testers with Memory Option

Features	SECUSTORE (Z745U)	SECUTEST SI (M702F)	SECUTEST SI+ (M702G)	SECUTEST PSI (GTMS016000R0001)	SECUTEST SIII+ Feature KB01 SECUIFEST	SECUTEST S2N + Option DBmed
Integrated printer for recording charts	—	—	—	•	—	—
Annotations via keyboard	—	•	•	•	—	—
Data memory (flash)	•	•	•	—	—	—
Data memory (battery buffered)	—	•	•	•	•	•
Protocol functions	•	•	•	•	—	—
Statistical evaluation of up to 8 instrument classes	—	•	•	•	—	—
Data transmission to PC via RS232 interface	•	•	•	•	•	•
Data transmission to PC via USB interface	—	—	•	—	—	—
Connection of a barcode scanner	•	•	•	•	•	•
Connection of an RFID scanner	•	•	•	•	•	•
Storage of function test values	•	•	•	•	—	—
Storage of comments on DUT	—	•	•	•	—	—

# SECUTEST SIII+ | . . .

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

### SECU-cal 10 Calibration Adapter

The calibration adapter is used for testing the measuring uncertainty of test instruments in accordance with DIN VDE 0701-0702 and IEC 62353 (VDE 0751-1). As a rule, these instruments must be tested once each year, as set forth by accident prevention regulation BGV A3 (previously VBG 4) as well as for certification in accordance with the ISO 9000 quality standard.




All limit values for the required tests per DIN VDE, as well as protective conductor resistance, insulation resistance, equivalent leakage current, differential and/or contact as well as housing leakage current must be tested.

### Test Adapter SECULOAD

Test Adapter for Testing Open-Circuit Voltage at Welding Units per EN 60974

In combination with a multimeter, the test adapter is used for testing welding units in accordance with the EN 60974-4:2007 standard. This standard stipulates that peak values for open-circuit voltage may not exceed the limit values regardless of the utilized settings.

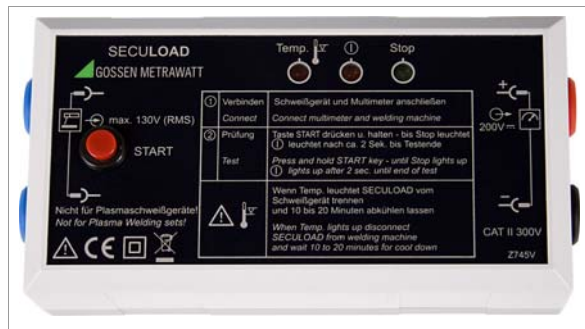

### K2010 Accessory Case for SECUTEST SIII+. . . and Accessories (not suited for Feature F02 or SECUTEST SIII+ H)



### F2000 Accessory Pouch for SECUTEST SIII+. . . and Accessories (not suited for Feature F02 or SECUTEST SIII+ H)





# Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

## Order Information

Designation	Type	Article Number
<b>Basic device</b>		
Test instrument with automatic test sequence, interface, German online instructions, earthing contact plug and socket, probe cable with test probe, plug-on alligator clip, 3 plug-on quick-connect terminals, DKD calibration certificate, operating instructions. See table on page 6 for features and expansions.	<b>SECUTEST SIII+ . . .</b>	<b>M7010</b> (all features: 00)
<b>Standard type available from stock</b>		
Test current ±200 mA DC or 25 A AC High-voltage test up to 6 kV DC Sequences for IEC 61010, IEC 60335, IEC 60950, data memory for up to 125 tests	<b>SECUTEST SIII+ H</b>	M7010-V013
<b>PC Analysis Software</b>		
For further information on software, please refer to our website  <a href="http://www.gossenmetrawatt.com">http://www.gossenmetrawatt.com</a> (→ Products → Electrical Testing → → Testing of Electric Appliances → SECUTEST ...)  or  <a href="http://www.gossenmetrawatt.com">http://www.gossenmetrawatt.com</a> (→ Products → Software → Software for Testers)		
<b>Accessories for Report Generation</b>		
SI module with RS232 and USB interface, with user languages D, GB, F, NL, I, E and CZ, batteries and operating instructions	SECUTEST SI+ <sup>D)</sup>	M702F
same as SECUTEST SI+, without USB interface, however, with additional integrated printer including 2 rolls of recording charts and 1 printer ribbon cartridge	SECUTEST PSI <sup>D)</sup>	GTM5016000R0001
Memory adapter for „direct print-out“ and internal test reports	SECUSTORE <sup>D)</sup>	Z745U
Firmware upgrade for SECUTEST data base: data memory for up to 125 tests (without function test values and without comments on DUT)	DBmed	Z853H
For barcode scanner, printer and RFID scanner see separate datasheet ID systems		
<b>Accessory Probes, Sensors, Adapters and Cables</b>		
Test probe with cable (no coil cord), 2 m, suitable for high-voltage test	SK2	Z745D
Test probe with cable (coil-cable), 2 meters long, suitable for high-voltage test	SK2W	Z745N
Probe cable 5 m	SK5	Z745K
Brush probe	Z745G	Z745G
12 conductor patient connection cable, each conductor with 4 mm plug for test instruments with feature J01	PA4	Z745L
Pt100 temperature sensor for surface and immersion measurements, -40 ... +600 °C	Z3409	GTZ3409000R0001
Pt100 oven sensor, -50 ... +550 °C	TF550	GTZ3408000R0001
Clip-on current sensor, can be set to 1 mA to 15 A or 1 A to 150 A, Frequency range: 45 ... 65 ... 500 Hz, 1 mV/mA and 1 mV/A	WZ12C <sup>D)</sup>	Z219C
Shunt for measuring range matching when using the instrument with feature G01 in combination with WZ12C transformer	Z864A	Z864A

Designation	Type	Article Number
Adapter for testing single-phase extension cables including earthing contact and inlet plug inserts, <b>do not use for high-voltage tests</b>	EL1	Z723A
Plug insert for EL1 in Switzerland per SEV	PRO-CH	GTZ3225000R0001
Plug insert for EL1 in Great Britain	PRO-GB	GTZ3226000R0001
Plug insert for EL1 GB measurement	PRO-GB/ring	GTZ3226000R0002
Plug insert for EL1 in Italy per IMQ	PRO-I	GTZ3227000R0001
Plug insert for EL1 in Denmark	PRO-DK	GTZ3219000R0001
Plug insert for EL1 in South Africa	PRO-RSA	Z501A
Plug insert for EL1 with 3 connector cables for any desired connection standards	PRO-UNI	GTZ3214000R0003
Plug insert for EL1 with 10 m cable for PE measurements and the like	PRO-RLO	GTZ3214000R0002
Plug insert Schuko or the like (replacement plug, included in EL1)	PRO-Schuko	GTZ3228000R0001
test adapter with single-phase and three-phase plug connectors up to CEE 32A – for all tests on single-phase and three-phase electrical devices without mains voltage per DIN VDE – for tests on single-phase and three-phase extension cables per DIN VDE	VL2E	Z745W
3 phase 16 A differential current adapter	AT16-DI	Z750A
3 phase 32 A differential current adapter	AT32-DI	Z750B
Test adapter for tests on devices with CEE16 and CEE32 connections (load rating of max 20 A)	AT3-II-S <sup>D)</sup>	Z745T
same as AT3-II-S, however, with a load rating of 32 A	AT3-II S32 <sup>D)</sup>	Z745X
3-phase current adapter 16A/32A (test case) for connection to the test instrument for tests per DIN VDE 0701-0702/IEC 62353 (VDE 0751) and IEC 601	AT3-III-E <sup>D)</sup>	Z745S
Test adapter in combination with <b>SECUTEST SIII+ . . .</b> for testing welding equipment in accordance with EN 60974-4:2007 standard.	SECULOAD	Z745V
Adapter for connecting devices under test: 3-pole 16 A, 5-pole 16 A and 32 A, 5 ea. 4 mm jack – for all tests in accordance with DIN VDE without line voltage at single and 3 phase electrical devices	CEE-Adapter	Z745A
Cable set for connecting test instruments to the mains without earthing contact socket and for connection of DUTs, <b>do not use for high-voltage tests</b>	KS13	GTY3624065P01
Cable set (1 pair of measuring cables) 1.2 m, with VDE-GS sign1000 V/CAT III, 600 V/CAT IV 16 A	KS17-2	GTY3520034P01
<b>Further Accessories</b>		
Calibration adapter for test instruments per DIN VDE 0701-0702 and IEC 62353 (VDE 0751) (max. 200 mA), <b>do not use for high-voltage tests and for protective conductor test current of 10 A or 25 A</b>	SECU-cal 10	Z715A
Pouch for SECUTEST... w/o HV module	F2000 <sup>D)</sup>	Z700D
Pouch big for SECUTEST with HV module and accessories	F2020	Z700F
Carrying case SECUTEST... w/o HV module	K2010	Z504L

<sup>D)</sup> Data sheet available

For additional information on accessories, please refer to:

- Measuring Instruments and Testers catalog
- our website [www.gossenmetrawatt.com](http://www.gossenmetrawatt.com)

# SECUTEST SIII+ | . . .

## Test Instrument for DIN EN 60601/60335/60950/61010, DIN VDE 0700/0701-0702 and IEC 62353 (VDE 0751-1)

---



Παγγαίου 6, Ν. Φιλοθέη Αμαρουσίου

Αθήνα 15123

Τηλ: 210 6754801, Fax: 210 6754804

[info@enia.gr](mailto:info@enia.gr)

[www.enia.gr](http://www.enia.gr)

---

Edited in Germany • Subject to change without notice • A pdf version is available on the internet

 **GOSSEN METRAWATT**

GMC-I Messtechnik GmbH  
Südwestpark 15  
90449 Nürnberg • Germany

Phone +49 911 8602-111  
Fax +49 911 8602-777  
E-Mail [info@gossenmetrawatt.com](mailto:info@gossenmetrawatt.com)  
[www.gossenmetrawatt.com](http://www.gossenmetrawatt.com)